

# HYDRIM® G4

C61W G4 Instrument Washer



## Automated Instrument Washing

the first step in the sterilization process



### Efficient

The time saved by automated washing compared to manual cleaning is approximately 1 hour of labor for every nine procedural set-ups<sup>1</sup>. HYDRIM® G4 cycle times are as fast as 40 minutes from start to dry. HYDRIM® ensures that every set of instruments receives the same validated cleaning process.



### G4 Technology

There is nothing like SciCan's G4 Technology. The large touch screen gives in-depth cycle information and allows users to view current operations.

From the networking portal, users can access cycle history and validation information. If there is ever a cycle fault, the G4 software will send an email to the office and the service technician, allowing remote troubleshooting.



### Capacity

HYDRIM® C61W G4 has a 6 full-size cassette capacity, allowing offices to turnover 120 instruments in as little as 40 minutes, start to dry.



<sup>1</sup> Tholen, Mark. "Increase production and decrease expenses in sterilization." Dental Economics, Volume 99, Issue 10, October 1, 2009.

## HYDR/M C61W G4 Specifications

### Unit Size

Dimensions (W x H x D):	598 x 520 x 526mm / 23.5" x 20.5" x 20.7"
Depth with Door Open:	82.9 cm / 32.6"
Weight:	44 kg / 97 lb

### Cycle Times\*

Rinse and Hold	5 min
Regular Wash + Dry	26 min + 10 min
Heavy Duty Wash + Dry	30 min + 10 min

### Water Consumption Per Cycle

Rinse and Hold	7.6 L / 1.7 Gal
Regular Wash	11.5 L / 3 Gal
Heavy Duty Wash	11.5 L / 3 Gal

### Electrical Connection

North America	208-240 V, 12 A, 60 Hz
Running Noise	65 dB

### Water and Drain Connections

Hot and Cold Water Tap Connection	3/4" Garden Hose, Male
Recommended Hot Water Temperature	60° C / 140° F
Cold Water Temperature	< 30° C / 86° F
Inlet Water Pressure	1-10 bar / 14.5-145 psi
Drain	2 cm / 3/4"

\*Actual cycle times will vary with incoming water temperature and pressure.

## Recommended Accessories

Part Number	Description
01-113251	6-Cassette Rack
01-113252	4 Cassette Rack
01-113253	3-Cassette, 2-Basket Rack
01-113254	10-Tray, 2-Basket Rack
01-113257	Trolley Center Support
01-113545	Large Basket
01-113546	Hygiene Basket
01-113547	Basket with Hinged Lid
CS-HIPC-U	HIP ULTRA Cleaning Solution



HYDR/M L110W G4  
also available

## VistaPure™ Water Filtration


Provides two grades of high-quality water for autoclaves, dental bottles, ultrasonic cleaners and the rinse cycle of instrument washers.

V3000

## HYDR/M VistaPure™ Connection

Protect expensive instruments and equipment by using the VistaPure system's hyperfiltered water for the final-rinse cycle of the HYDR/M for a spot-free finish.

V3010

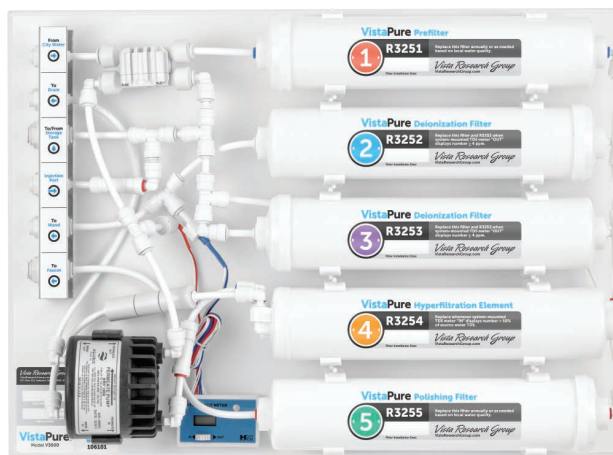

VistaPure™ Water Filtration shown above.

Manufactured by:  
**SciCan Ltd.**  
1440 Don Mills Rd.,  
Toronto, ON, Canada, M3B 3P9  
T +1 416-445-1600  
TF +1 800-667-7733  
customerservice@scican.com

Distributed by:  
**SciCan Inc.**  
701 Technology Dr.,  
Canonsburg, PA 15317, USA  
T +1 724-820-1600  
TF +1 800-572-1211  
uscustomerservice@scican.com

VistaPure™, VistaCool™ and DrainHub™  
are manufactured by:  
**Crosstex International, Inc.**  
1733 W. Parkside Lane  
Phoenix, AZ 85027 / USA  
T +1 800-531-2221  
VistaResearchGroup.com

**COLTENE**